



## PRODUCT CATALOGUE

Soft Lying Hyperbaric Chamber

Soft Sitting Hyperbaric Chamber

Wheelchair Hyperbaric Chamber

Multiplace Hyperbaric Chamber

Hard Type Hyperbaric Chamber

2.0ATA Hard Hyperbaric Chamber

Stretcher Type Hyperbaric Chamber

Veterinary Hyperbaric Chamber



## COMPANY INTRODUCTION



Shanghai Baobang Medical Equipment Co.,Ltd, located in Shanghai, China, is a professional manufacturer, specialized in the development, production and export of hyperbaric chambers. The company works under ISO13485, that represents the requirements for a comprehensive quality management system for the design and manufacture of this product. Presently the company has succeeded in the production of serial hyperbaric chambers from soft and hard lying type, soft sitting type, to veterinary hard type chambers, which have been widely applied in homecare, clinics, wellness centers, beauty salons, spas, sports industry, and high altitude area, etc. In the past 17 years, MACY-PAN hyperbaric chambers are sold over 126 countries and regions, enjoy a high reputation in international market.

We present the strong quality commitment and we have gone through hard work to be able to offer this health equipment at an accessible price. We are proud to be the top Chinese factory doing this chamber and be able to give an opportunity to people that need HBOT in their lives.

The basic business policy of the company has always meant a very high standard, reliable and high quality service. We welcome not only family buyers, but also sales companies by providing you with tailored OEM products. With confident and enthusiastic service, we expect our cooperation with products of sound quality and affordable price.





## CERTIFICATE



**INGERCERTIFICATION ASSESSMENT SERVICES**

**MANAGEMENT SYSTEM CERTIFICATE**

Certificate No.: 117 17 2003 004 423  
This is to certify the quality management systems of

**Shanghai Baobang Medical Equipment Co., Ltd.**

Unified Social Credit Code: 91310117667825979A

Location: Building 15, No. 889 Guilan Road, Shishuang Town, Songjiang District, Shanghai

has been assessed and registered as meeting the requirements of GB/T 19001-2016/ISO 9001:2015

Scope of approval:  
Design, Production and Export of Portable Pressure Chamber, Hand Structure Pressure Chamber and Sauna Cabin

Signed by: [Signature]

Approved Date: 10 Dec. 2017  
Expiry Date: 10 Feb. 2021

Shanghai Inger Certification Assessment Co., Ltd.  
100000 Shanghai, China

**INGERCERTIFICATION ASSESSMENT SERVICES**

**MANAGEMENT SYSTEM CERTIFICATE**

Certificate No.: 117 17 2003 004 423  
This is to certify the quality management systems of

**Shanghai Baobang Medical Equipment Co., Ltd.**

Unified Social Credit Code: 91310117667825979A

Location: Building 15, No. 889 Guilan Road, Shishuang Town, Songjiang District, Shanghai

has been assessed and registered as meeting the requirements of GB/T 19001-2016/ISO 9001:2015

Scope of approval:  
Design, Production and Export of Portable Pressure Chamber, Hand Structure Pressure Chamber and Sauna Cabin

Signed by: [Signature]

Approved Date: 10 Dec. 2017  
Expiry Date: 10 Feb. 2021

Shanghai Inger Certification Assessment Co., Ltd.  
100000 Shanghai, China

**Certificate of Compliance**

No. 2017112518M929E

Certificate's Holder: Shanghai Baobang Medical Equipment Co., Ltd.  
Building 15, No. 889 Guilan Road, Shishuang Town, Songjiang District, Shanghai

Certification ECM Mark:

Product: Oil-free Oil Compressor  
Model(s): MP92524, MP92524.1, MP92524, BK101, K101, BK106, K106, BDF03, BDF03, BDF100

Verification to: Standard: EN ISO 12100:2010  
EN 60204-1:2006+A1:2009+A12:2010, EN 61010-1:2010/AC:2012  
Related to CE Directives: 2004/AC/EC (Machinery), 2014/53/EU (Low Voltage), 2014/53/EU (Electromagnetic Compatibility)

Remarks: The product(s) has been verified on a voluntary basis. The product(s) satisfies the requirements of the Certification mark of CE mark in the above listed standards. The above Compliance mark can be placed on the product(s) according to the CE marking rules to indicate that the product(s) has been assessed and registered as meeting the requirements of the standards listed in the scope of approval. The Certificate of Compliance can be placed on the product(s) to indicate that the product(s) has been assessed and registered as meeting the requirements of the standards listed in the scope of approval.

Issue Date: 10 Dec. 2017  
Expiry Date: 10 Feb. 2021

Shanghai Inger Certification Assessment Co., Ltd.  
100000 Shanghai, China

**Certificate of Compliance**

No. 2017122518M929E

Certificate's Holder: Shanghai Baobang Medical Equipment Co., Ltd.  
Building 15, No. 889 Guilan Road, Shishuang Town, Songjiang District, Shanghai

Certification ECM Mark:

Product: Portable Bag  
Model(s): ST701, ST602, ST603, ST702, ST703, ST704, ST705, ST706, ST707, ST708, ST709, ST710, ST711, ST712, ST713, ST714, ST715, ST716, ST717, ST718, ST719, ST720

Verification to: Standard: EN 13557:2002

Related to CE Directives: 2001/95/EC (General Product Safety)

Remarks: The product(s) has been verified on a voluntary basis. The product(s) satisfies the requirements of the Certification mark of CE mark in the above listed standards. The above Compliance mark can be placed on the product(s) according to the CE marking rules to indicate that the product(s) has been assessed and registered as meeting the requirements of the standards listed in the scope of approval. The Certificate of Compliance can be placed on the product(s) to indicate that the product(s) has been assessed and registered as meeting the requirements of the standards listed in the scope of approval.

Issue Date: 10 Dec. 2017  
Expiry Date: 10 Feb. 2021

Shanghai Inger Certification Assessment Co., Ltd.  
100000 Shanghai, China

**Certificate of Compliance**

No. 2018011320M929E

Certificate's Holder: Shanghai Baobang Medical Equipment Co., Ltd.  
Building 15, No. 889 Guilan Road, Shishuang Town, Songjiang District, Shanghai

Certification ECM Mark:

Product: Stainless Steel Tank  
Model(s): H9700, H9710

Verification to: Standard: EN 1447:2002

Related to CE Directives: 2014/68/EU (Pressure Equipment)

Remarks: The product(s) has been verified on a voluntary basis. The product(s) satisfies the requirements of the Certification mark of CE mark in the above listed standards. The above Compliance mark can be placed on the product(s) according to the CE marking rules to indicate that the product(s) has been assessed and registered as meeting the requirements of the standards listed in the scope of approval. The Certificate of Compliance can be placed on the product(s) to indicate that the product(s) has been assessed and registered as meeting the requirements of the standards listed in the scope of approval.

Issue Date: 10 Dec. 2017  
Expiry Date: 10 Feb. 2021

Shanghai Inger Certification Assessment Co., Ltd.  
100000 Shanghai, China

**Verification of Compliance**

No. EC: 1202.00140515.58M0747  
National Commission No. EC-1202/001

Certificate's Holder: Shanghai Baobang Medical Equipment Co., Ltd.  
No. 796, Yuyang Road, Songjiang District, Shanghai, China

Product: Portable Hypertension Oxygen Chamber  
Model(s): ST701, SH602, ST701-3, ST702, ST703, ST704, ST705, ST706, ST707

Verification to: Standard: EN 14831:2009

Related to CE Directives: 2007/02/EC (Medical Devices)

Remarks: The Verification of Compliance has been issued if a holding body, ECN number 0014, has received the Declaration of Conformity (DoC) for the product(s) listed in the scope of approval. The Declaration of Conformity (DoC) is a statement by the manufacturer that the product(s) meets the requirements of the standards listed in the scope of approval. The Declaration of Conformity (DoC) is a statement by the manufacturer that the product(s) meets the requirements of the standards listed in the scope of approval. The Declaration of Compliance can be placed on the product(s) to indicate that the product(s) has been assessed and registered as meeting the requirements of the standards listed in the scope of approval.

Issue Date: 10 Dec. 2017  
Expiry Date: 10 Feb. 2021

Shanghai Inger Certification Assessment Co., Ltd.  
100000 Shanghai, China

外观设计专利证书

专利号: 2017 2 0588663.5

外观设计名称: 负压空气呼吸机 (结构)

设计人: 周志南

专利号: 2017 2 0588663.5

申请日期: 2017 年 11 月 27 日

专利权人: 上海宝邦医疗器械有限公司

授权公告日: 2018 年 04 月 20 日

周志南 申志南

外观设计专利证书

专利号: 2017 2 0588663.5

外观设计名称: 负压空气呼吸机 (结构)

设计人: 周志南

专利号: 2017 2 0588663.5

申请日期: 2017 年 11 月 27 日

专利权人: 上海宝邦医疗器械有限公司

授权公告日: 2018 年 04 月 20 日

周志南 申志南

发明专利证书

专利号: 2017 1 0588663.5

发明名称: 负压空气呼吸机 (结构)

发明人: 周志南

专利号: 2017 1 0588663.5

申请日期: 2017 年 11 月 27 日

专利权人: 上海宝邦医疗器械有限公司

授权公告日: 2018 年 04 月 20 日

周志南 申志南


# Hyperbaric Oxygen Therapy




Hyperbaric Oxygen Therapy is the medical use of oxygen at an increased level of pressure, (hyper) means increased, (baric) means pressure. Hyperbaric Oxygen Therapy is a treatment which enhances the body's natural healing process by providing an environment which allows the body to absorb much higher amounts of oxygen than possible at normal atmospheric pressure. It increases the oxygen transport capacity in the blood. This noninvasive therapy is the most trusted way to increase oxygen levels to all organs of the body.

Henry's Law:

**“The amount of dissolved gas in a liquid is proportional to its partial pressure above the liquid.”**

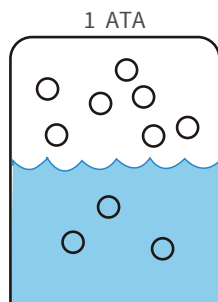


Circle: Gas  
Blue: Liquid

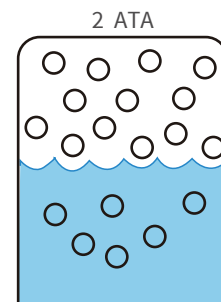

the state of gas dissolved in liquid at one ATA pressure.

the state of gas dissolved in liquid at two ATA pressure.

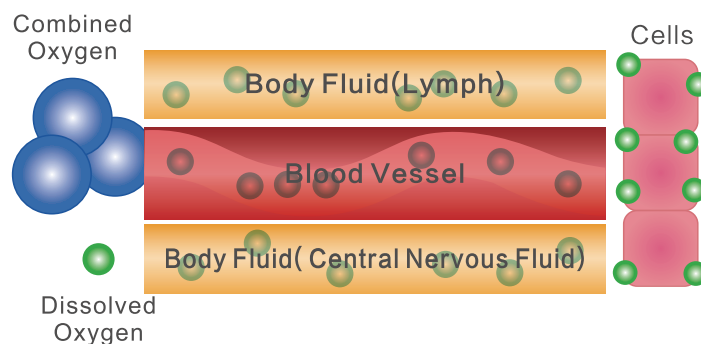

Oxygen transport can be classified into “combined oxygen and dissolved oxygen”. When entering blood, most oxygen will be combined with hemoglobin, then become combined oxygen. Another part of oxygen will be dissolved in plasma directly and physically, then become dissolved oxygen. Therefore, combined oxygen is difficult to pass the narrow vascular wall, while dissolved oxygen is small and can greatly dissolved in blood and fluid under hyperbaric circumstances.




### Clinical

For Medical: Clinic, Hospital, High Altitude Area  
 Reduces inflammation.  
 Speed Wound Healing.  
 Improves immune system support.  
 ALS, PTSD, Diabetes, Autism, Epilepsy treatment etc.



### Sports & Recovery

Eliminate sport fatigue  
 Regather strength  
 Alleviate sport damage  
 Rapidly remove lactic acid  
 Relieve damage caused by free radical

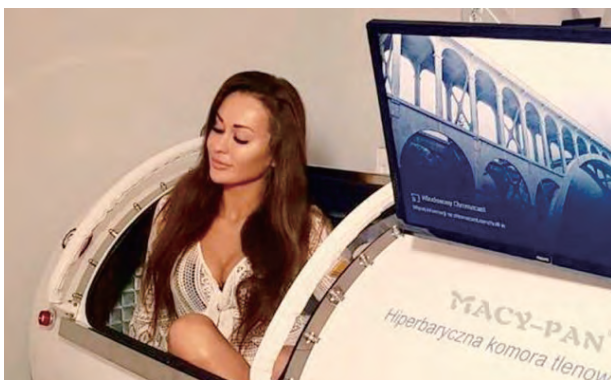

### Beauty&Anti-aging

Stimulate the natural recovery process and immunity.  
 Boost vitality, heighten concentration and memory.  
 Improve cell rejuvenation and help with anti-aging strategies.  
 Remove the sensation of brain fog and jet lag.  
 Eliminates toxins and debris.  
 Reactivate skin cells and improve blood circulation.



### Home care

Oxygen to stimulate the natural healing process  
 Hyperbaric oxygen therapy sessions are comfortable and relaxing and all you have to do is breathe oxygen inside a chamber with increased atmospheric pressure that enhances the distribution of oxygen through blood capillaries, through tissues, organs, and the brain. The treatment is non-invasive and fit for mostly everyone.

# Hyperbaric Chamber Lying Type



## Model:ST801 / ST901



### Features

- 1.3 ATA /1.4 ATA /1.5 ATA operating pressure.
- 2 zipper seal for easier entrance.
- 7 large transparent viewing windows to prevent claustrophobia.
- Cotton chamber protection cover, to avoid dirty and easy to wash.
- Most popular model for home treatment or commercial usage.
- Emergency valve - To speed the depressurization in emergency.
- Easy Operation - One people can operate it without assistance.
- Delivers 93% oxygen under pressure via an oxygen headset/facial mask.
- Internal Pressure Gauge - User can observe the pressure while inside the chamber.
- Internal metal frame - Hold the shape when deflated.

Model		ST801	ST901
Chamber Capsule	Pressure	1.3ATA /1.4ATA/1.5ATA	1.3ATA/1.4ATA
	Material	TPU	TPU
	Size ( D*L )	80x225cm (32x89 inch)	90x225cm (36x89 inch)
	Weight	13kg	16kg
Air Compressor	Size	39x24x26cm (16x10x11 inch)	
	Type	Quite&Oil Free Type.Two Super adsorption activated carbon filters in air Input and Outlet to protect.	
Oxygen Concentrator	Size	35x32x57cm (14x13x23 inch)	
	Flow Rate	5Liter/min /10Liter/min	
	Description	PSA Molecular Sieve High Technology.Continuouss Oxygen production,no need for oxygen tank.	
Air Cooler	Size	18x12x35cm(7.2x5x14 inch)	
	Description	Cooling and dehumidification	

## Hyperbaric Chamber Sitting Type



### Model:ST1700/ST2200



#### Features

1.3 ATA /1.4 ATA /1.5 ATA operating pressure.

Fully covered with nylon cover - Unique 3 zipper seal for air leak prevention.

Designed to add a folding chair(adjustable height) for ultimate comfort.

Internal metal frame - Hold the shape when deflated.

Emergency valve - To speed the depressurization in emergency.

Delivers 93% oxygen under pressure via an oxygen headset/ facial mask.

Easy Operation - One people can operate it without assistance.

Internal Pressure Gauge - The user can observe the pressure while inside the chamber.

Extended site sitting chamber,one can both lay down and sit down for ST2200.

Model		ST1700	ST2200
Chamber Capsule	Pressure	1.3ATA /1.4ATA/1.5ATA	1.3ATA/1.4ATA
	Material	TPU	TPU
	Size(D*L)	170x70x110cm (67x28x43 inch)	220x70x110cm(89x28x43 inch)
	Weight	13kg	14 kg
	Type	Sitting/Semi-Lying	Sitting/Lying
Air Compressor	Size	39x24x26cm (16x10x11 inch)	39x24x26cm (16x10x11 inch)
	Type	Quiet&Oil Free Type. Two Super adsorption activated carbon filters in air Input and Outlet to protect.	Quiet&Oil Free Type. Two Super adsorption activated carbon filters in air Input and Outlet to protect.
Oxygen Concentrator	Size	35x32x57cm (14x13x23 inch)	35x32x27cm(14x13x23 inch)
	Flow Rate	5Liter/min /10Liter/min	5Liter/min/10Liter/min
	Description	PSA Molecular Sieve High Technology. Continuous Oxygen production, no need for oxygen tank.	PSA Molecular Sieve High Technology. Continuous Oxygen production, no need for oxygen tank.
Air Cooler	Size	18x12x35cm(7.2x5x14 inch)	18x12x35cm(7.2x5x14 inch)
	Description	Cooling and dehumidification	Cooling and dehumidification

# Hyperbaric Chamber Sitting Type



## Model:L1



### Features

1.3 ATA /1.4 ATA /1.5 ATA operating pressure.

“L” shaped zipper for easier entrance.

Vertical design - small footprint.

Most economical model for home treatment or commercial usage.

Emergency valve - To speed the depressurization in emergency.

Easy Operation - One people can operate it without assistance.

Delivers 93% oxygen under pressure via an oxygen headset/facial mask.

Internal Pressure Gauge - User can observe the pressure while inside the chamber.

Internal metal frame - Hold the shape when deflated.

Model		L1
Chamber Capsule	Pressure	1.3/1.4/1.5 ATA
	Material	TPU
	Size(D*L)	140x100x160cm (56x40x64 inch)
	Weight	13kg
Air Compressor	Size	39x24x26cm (16x10x11 inch)
	Type	Quiet&Oil Free Type. Two Super adsorption activated carbon filters in air Input and Outlet to protect.
Oxygen Concentrator	Size	35x32x57cm (14x13x23 inch)
	Flow Rate	5 Liter/min /10Liter/min
	Description	PSA Molecular Sieve High Technology. Continuous Oxygen production, no need for oxygen tank.
Air Cooler	Size	18x12x35cm (7.2x5x14 inch)
	Description	Cooling and dehumidification



## Hyperbaric Chamber Wheelchair Type



### Model:MC4000



#### Features

- 1.3ATA mild operating pressure.
- Fully covered with nylon cover - Unique 3 zipper seal for air leak prevention.
- Emergency valve - To speed the depressurization in emergency.
- Delivers 93% oxygen under pressure via an oxygen headset or oxygen mask.
- Available "u" type door and "n" type for your choice
- Designed to accommodate 2 folding chairs with comfortable space
- Wheelchair accessible, Designed for people with limited mobility.

Model		MC4000
Chamber Capsule	Pressure	1.3/1.4 ATA
	Material	TPU
	Size(D*L)	140x130x175cm (56x52x70 inch)
	Weight	20 kg
Air Compressor	Size	39x24x26cm(16x10x11 inch)
	Type	Quiet&Oil Free Type. Two Super adsorption activated carbon filters in air Input and Outlet to protect.
Oxygen Concentrator	Size	35x32x57cm(14x13x23 inch)
	Flow Rate	5 Liter/min /10Liter/min
	Description	PSA Molecular Sieve High Technology. Continuous Oxygen production, no need for oxygen tank.
Air Cooler	Size	18x12x35cm(7.2x5x14 inch)
	Description	Cooling and dehumidification
Folding Chair (Optional)	Description	Durable and elastic. Can sit or semi-lying inside chamber. Weight: 9kg

# Hyperbaric Chamber Multiplace Type



## Model:STM2000



### Features

1.3 ATA operating pressure.

Diameter is 2000mm. Could accommodate 4-6 people.

Special material design at the bottom, smooth and flat. Could sit inside or lie down

Emergency valve - To speed the depressurization in emergency.

Delivers 93% oxygen under pressure via an oxygen headset/ facial mask.

Easy Operation - One people can operate it without assistance.


Internal Pressure Gauge - The user can observe the pressure while inside the chamber.

Model		STM2000
Chamber Capsule	Pressure	1.3ATA
	Material	TPU
	Size(D*L)	200x180x200cm(78x70x78 inch)
	Weight	60 kg
	Type	Multiplace
Air Compressor	Size	76x42x66cm(30x16x26 inch)
	Type	Quiet&Oil Free Type. Two Super adsorption activated carbon filters in air Input and Outlet to protect.
Oxygen Concentrator	Size	35x32x27cm(14x13x23 inch)
	Flow Rate	10 Liter/min
	Description	PSA Molecular Sieve High Technology. Continuous Oxygen production, no need for oxygen tank.
Folding Chair (4pcs)	Description	Durable and elastic. Can sit or semi-lying inside chamber. Weight: 9kg

## Hyperbaric Chamber Multiplace Type



### Model: HE5000

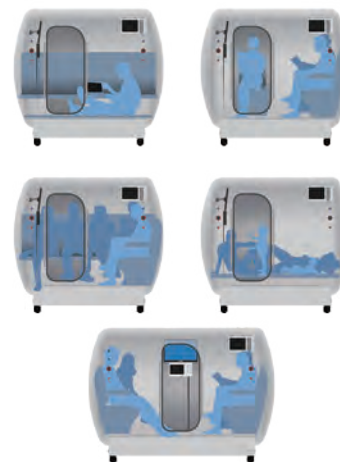

#### Features

**1.5/2 ATA operating pressure**

**Emergency valve - To speed the depressurization in emergency.**

**Delivers 93% oxygen under pressure via an oxygen headset or oxygen mask.**

**Internal Pressure Gauge - The user can observe the pressure while inside the chamber.**

**Specially designed for 3-6 people use, people can sit/semi-sitting inside chamber with friends or families and enjoy hyperbaric oxygen therapy together.**

**More space for people who are suffering from claustrophobic.**

**The bracket at the bottom can be used to stabilize the chamber.**

Model		HE5000
Chamber Capsule	Pressure	1.5ATA/2ATA
	Material	304 stainless steel
	Size(D*L)	210x180x176cm (83x71x69 inch)
	Weight	480 kg
System Machine	Size	172x35x65cm (38x13x24 inch)
	Flow of Air	120L/Minute
	Power	1000W
	Voltage	110/220/240V
Oxygen Concentrator *2	Size	47x37x59cm (19x15x24 inch)
	Flow	20L/min (adjustable)
	Power	500W
	Voltage	110/220/240V

# Hyperbaric Chamber Hard Type



## Model:HP1501



### Features

1.5/1.6 ATA operating pressure

75 and 90cm diameter optional.

Interphone system for two-way communications.

Pneumatic control system, unique sliding door secure locking mechanism for convenient door opening and closing.

Wide and convenient transparent viewing window.

Delivers 93% oxygen under pressure via an oxygen headset/facial mask.

User friendly design, safe and simple to operate.

Control system combines of air compressor, oxygen concentrator and cooling unit.

Model		HP1501-75	HP1501-90	HP1501-100
Chamber Cabin	Pressure	1.5ATA/1.6ATA	1.5ATA/1.6ATA	1.5ATA/1.6ATA
	Material	Stainless steel+Polycarbonate		
	Size ( D*L )	75x220cm (30x87 inch)	90x220cm (36x87 inch)	100x220cm (39x87 inch)
	Weight	140kg	170kg	220kg
System Machine	Size/Weight	77x46x72cm (30x17x28 inch)		
	Flow of Air	72 Liter/min		
	Flow of Oxygen	5Liter/min /10Liter/min		
	Air Filter	Dual		
	Voltage	110VAC/220VAC ± 10%		
	Description	All-in-one machine, a combine of air compressor, oxygen concentrator, and cooler		



## Hyperbaric Chamber 2.0ATA Type



### Model:HP2202



#### Features

2.0 ATA operating pressure.

The highest pressure in MACY-PAN chambers.

Pneumatic control system, unique sliding door secure locking mechanism for convenient door opening and closing.

75 and 85cm diameter optional.

Interphone system for two-way communications.

Delivers 93% oxygen under pressure via an oxygen headset/facial mask.

Control system combines of air compressor oxygen concentrator and cooling unit.

User friendly design, safe and simple to operate.

Model		HP2202-75	HP2202-85	HP2202-90
Chamber Cabin	Pressure	2.0ATA	2.0ATA	2.0ATA
	Material	Stainless steel+Polycarbonate		
	Size ( D*L )	75x220cm (30x87 inch)	85x220cm (34x87 inch)	90x220cm (36x87 inch)
	Weight	150kg	180kg	200kg
System Machine	Size/Weight	76x45x65cm (30.4x18x26 inch), 80kg		
	Flow of Air	72 Liter/min		
	Flow of Oxygen	5Liter/min / 10Liter/min		
	Air Filter	Dual		
	Voltage	110VAC/220VAC ± 10%		
	Description	All-in-one machine, a combine of air compressor, oxygen concentrator, and air dryer, air conditioner as an option		

# Hyperbaric Chamber Stretcher Type



## Model:MC750



### Features

1.8ATA operating pressure.

Patient with limited mobility can be moved in and out of the chamber via a sliding stretcher.

Delivers 93% oxygen under pressure via an oxygen headset or oxygen mask.

Communication system allows the patient to communicate to the operator.

Dual-operated deflate valve, allows the professional to regulate them from the outside, or the patient to operate them internally.


The Control system combines compressor, oxygen concentrator, and cooling unit, pressure, temperature, light within the chamber is easily controlled by the operator.

Model		MC750
Chamber Cabin	Pressure	1.8ATA
	Material	Stainless steel
	Size ( D*L )	75x220 cm (30x87 inch)
	Weight	132kg
System Machine	Size/Weight	58x42x56cm (23x17x22 inch)
	Flow of Air	80L/min
	Flow of Oxygen	5Liter/min /10Liter/min
	Air Filter	Dual
	Voltage	110/220/240VAC ± 10%
	Description	A combine of air compressor, oxygen concentrator, and air cooler
Air Conditioner	Temperature range	110°C~Room temperature
	Voltage	110/220/240VAC ± 10%
	Weight	15kg

## Hyperbaric Chamber Veterinary Type



### Model:HP600 / HP700

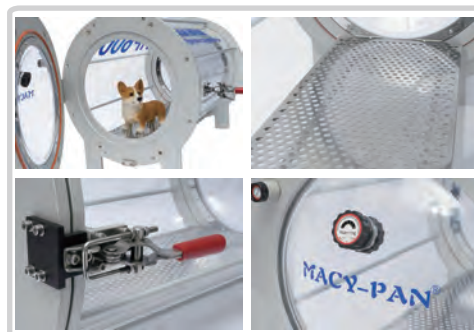

#### Application Fields

Gastric dilation and volvulus, and the associated injuries.

Fracture repair, especially open fractures or those involving malunion or nonunion

Gastrointestinal resection and anastomosis, especially post-foreign body obstruction.

Amputations

Skin grafts

Hemilaminectomy

Ear canal ablation

Model		HP600	HP700
Chamber Cabin	Pressure	1.3/1.4ATA(4/5,8psi)	1.5ATA(7psi)
	Material	Aluminum+Acrylic	Stainless Steel +PC
	Cabin diameter	60cm	70cm
	Size(D*L)	125x73x86cm (49x29x34 inch)	142x71x93cm (57x28x37 inch)
	Weight	65kg	100kg
System Machine	Size/Weight	77x46x72cm (30x18x28 inch),70kg	
	Flow of Air	72 Liter/min	
	Flow of Oxygen	5Liter/min /10Liter/min	
	Air Filter	Dual	
	Voltage	110VAC/220VAC ± 10%	
	Description	All-in-one machine, a combine of air compressor, oxygen concentrator, and cooler	



Sold around the world







### Why is oxygen so important?

Every day an average adult consumes 2 kg of food, 1 liter of water and almost 550 liters of oxygen. 90 % of our daily energy comes from oxygen and the rest food and water. Human brain is 2 % of the total weight and uses 25% of the oxygen we breath. We need this oxygen for the energy cycle that sustains life. When we do not have enough oxygen in our body tissues a series of events occur that if not corrected lead to disease conditions, either infection, tissue destruction or both. If there is low oxygen in tissues (hypoxia) there is a short window of opportunity to correct it. An excellent method to correct tissue hypoxia is by using a hyperbaric chamber.

### What do you feel inside a hyperbaric chamber?

When the chamber is pressurization, your ears may feel pressure changes. You might feel a little painful in ears. To equalize the pressure and avoid the feeling of fullness in ears, You can yanwing, swallowing or "blow the nose". Other than this ear pressure there are no unusual or different sensations.

### How does being inside a pressurized chamber give us more oxygen?

When we are inside a chamber pressurized at 1.3 or 1.5 times of the normal air pressure it may not feel different, but breathing pure oxygen in such a chamber gives us 7-8 times the regular amount of oxygen. Henry's Law- Henry's law states that for a gas-liquid interface the amount of the gas that dissolves in the liquid is proportional to its partial pressure. The pressurized oxygen dissolves directly into all body fluids. In a few minutes this extra oxygen builds up tissue oxygen levels far above normal. This action has been scientifically proven to stimulate healing.

### What improvements can I expect to see ?

This is one of those questions which has so many different variables and answers we will have to generalize. We will answer this question from general experience and with the understanding that no 2 people get the same results. It is very rare that we do not see improvements. What we are looking for is an improvement in the quality of life for anyone.

### When will I know how many sessions of hbot are enough?

You will not. What you will be able to notice is whether you are continuing to get improvements. A clinical study on 1730 subjects shows the need of around 15 sessions for acute injuries and 40 for chronic ones. If you saw results during the first 40 sessions, do another 40. If improvements continue, do another 40. Etc. Etc. You are a better judge of the number of sessions needed than any Doctor.

### What about Oxygen Toxicity?

Oxygen Toxicity is not a problem at pressures below 1.5 ATA. Many times stress is mistaken for oxygen toxicity. Occasionally, there may be a patient who is more sensitive to Oxygen than others and pressures may need to be adjusted to meet the needs of the patient.



**BETTER HEALTH THROUGH HBOT!**

**MACY-PAN<sup>®</sup>**

---

Shanghai Baobang Medical Equipment Co., Ltd  
Add: Bldg.15, No.889, Guinan Rd., Songjiang Dc., Shanghai, China 201617.  
E-mail: [rank@macy-pan.com](mailto:rank@macy-pan.com)  
Tel: 0086 18017021619  
Web: [www.hbotmacypan.com](http://www.hbotmacypan.com)